

medisana®

MG 200

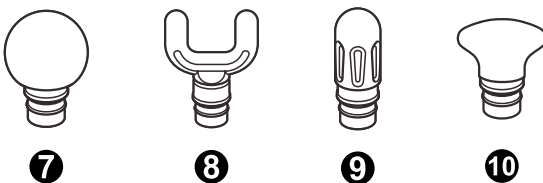Manual

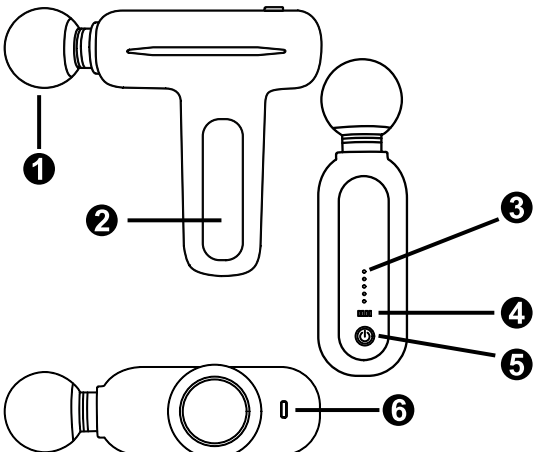Service

medisana GmbH
Carl-Schurz-Str. 2,
41460 NEUSS, Germany

88998

40 15588 88998 1

Handmassagegerät MG 200

DE

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Handheld Massager MG 200

GB


- medisana MG 200** (with holder, so that the massage head can be easily exchanged)
- Handle
- LEDs displaying the vibration intensity
- LEDs displaying operation and battery charging level
- On/off button, and change vibration intensity
- USB-C port (charging socket)
- Round head – suitable for large muscle groups
- U-shaped head – suitable for the muscles in the area of the spine
- Cone-shaped head – suitable for massaging deep tissue
- Flat head – suitable for relaxation and beauty massage

Safety Instructions

Read the instruction manual carefully before using this device, especially the safety instructions, and keep the instruction manual for future use. These brief instructions are an integral part of the instructions for use. The complete instructions for use are available for download at <https://www.medisana.com>. Should you give this device to another person, it is vital that you also pass on the instructions for use.

- Please ensure that the polythene packing is kept away from the reach of children! Risk of suffocation!
- Only use the appliance for the purpose for which it is intended in accordance with the instructions for use. Any misuse will invalidate the warranty.
- The unit is not intended for commercial purposes or the medical sector.
- Only charge the unit's battery using the USB charging cable supplied. Never leave the unit unattended when it is connected to the mains (for charging) or in operation.
- Never charge the device in a flammable environment (e.g. flammable objects, gasoline, flammable gases, paint, etc.)
- Only connect the unit for charging so that the cable and power connector are freely accessible.
- The built-in batteries are not intended for replacement or removal by the end user. If problems occur with the power supply, contact the customer service.
- This device must only be supplied at safety extra low voltage (SELV) corresponding to the marking on the appliance.
- To disconnect the unit from the mains, always pull the plug out of the power connection. Do not pull on the cable!
- If you have any health concerns, consult your doctor before use.
- You must not use the appliance if you suffer from one or more of the following diseases or medical conditions: Circulatory disorders, varicose veins, aneurysms, open wounds, heart disease, bruises, skin tears, cancer, phlebitis or thrombosis.
- Do not use the device if you are pregnant. Special care is necessary when the appliance is used on, by or near children, sick or helpless persons.
- Always consult your doctor before use if you are taking blood thinners, wear a pacemaker or defibrillator, or have undergone internal surgery to assess possible risks of using the device.
- If you experience pain or find the massage uncomfortable, discontinue use and talk to your doctor.
- In case of unexplained pain, if you are undergoing medical treatment and/or using medical equipment, consult your doctor before using the massager.
- Do not use this device to support or replace medical applications. Chronic conditions and symptoms could be aggravated.
- Keep the device out of the reach of children under 16 years.

- This appliance can be used by children aged 8 years and above and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge if they have been given supervision or instruction concerning use of the appliance in a safe way and understand the hazards involved.
- Children must not play with the appliance.
- Check the massager carefully for damage before each use. Do not operate a defective appliance.
- Do not use the appliance if any damage is visible on the appliance, if it is not working properly, if it has been dropped or if it has become wet. To avoid hazards, send the unit to the service centre for repair.
- Do not use the unit when you are sleeping or in bed.
- Do not use the appliance before falling asleep. The massage has a stimulating effect.
- Never use the appliance when driving a vehicle or operating machinery.
- Too long / too frequent use can cause overheating. Always allow the appliance to cool down before using it again.
- Never place or use the appliance directly next to an electric heater or other heat sources.
- Do not use the appliance on sensitive areas of the body, such as the head, face, spine, shins, backs of feet or hands, etc.
- This appliance should not massage the skin surface directly to avoid damage to the skin. The massage part should be covered with clean, dry clothes or a clean, dry towel. It should be pressed and moved gently. The massage time on the same spot should not exceed 60 seconds.
- Avoid contact of the appliance with pointed or sharp objects.
- Do not cover the appliance when it is switched on. Do not use the appliance under blankets or pillows. There is a risk of fire, electric shock and injury.
- Make sure that fingers or other parts of the body never come between the massage head holder and the massage heads or directly onto the massage head when switching on. Risk of injury!**
- Do not immerse the appliance in water or other liquids.
- However, if liquid enters the appliance, switch it off immediately and disconnect the power supply if it was connected. Wait until the unit is completely dry.
- Always switch off the appliance after use.
- Never block the ventilation slots of the device.
- Do not change the massage head while the device is switched on.
- Do not use the device on these parts of the body: head, groin area, chest, and do not use it on parts that have undergone a surgery.
- The unit is maintenance-free. You may only carry out cleaning work on the unit yourself. In the event of malfunctions or damage, do not repair the unit yourself, as this will invalidate any warranty claims. Consult your specialist dealer and have repairs carried out only by authorised service centres in order to avoid hazards.
- Only use accessories intended for the unit and supplied by the manufacturer.
- Cleaning and user maintenance must not be carried out by children without supervision.
- Voltage supply:** Li-Ion battery 7.4V, 1800mAh. Input: 5V== 2.0A
Charging time: approx. 3.5 hours

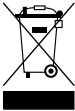Disposal: This product must not be disposed of together with domestic waste. All users are obliged to hand in all electrical or electronic devices, regardless of whether or not they contain toxic substances, at a municipal or commercial collection point so that they can be disposed of in an environmentally acceptable manner. Please remove the batteries before disposing of the device / unit. Do not dispose of old batteries with your household waste, but at a battery collection station at a recycling site or in a shop. Consult your local authority or your supplier for information about disposal.

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